



# Serani™ Max

## Analyzer Interface Software

GASES & CHEMICALS

CEMS

ENERGY

SEMI & HB LED

ATMOSPHERIC


LAB & LIFE SCIENCE

### Maximize Convenience!

Tiger Optics' Serani™ Max Interface Software offers convenient remote access from your computer or laptop for Tiger's latest generation analyzers, such as the HALO KA Max and T-I Max. The Windows-based software package communicates with the analyzer through Ethernet and allows full control of the system conveniently from your desk. It also offers a simple way to download and save measurement data and settings to your computer without requiring direct access to the analyzer.


### Serani™ Max features in detail:

- Monitoring and control of your Tiger analyzer remotely with dedicated application software installed on your computer or laptop
- Data plotting and download to your computer
- Easy change of analyzer settings, such as gas type and data filter
- One-step data collection for "Performance Verification"
- Compatible with Windows 7, 8 and 10
- Compatible with any "Max" series analyzer, as well as select HALO, ALOHA H<sub>2</sub>O and Spark series analyzers (details see reverse)

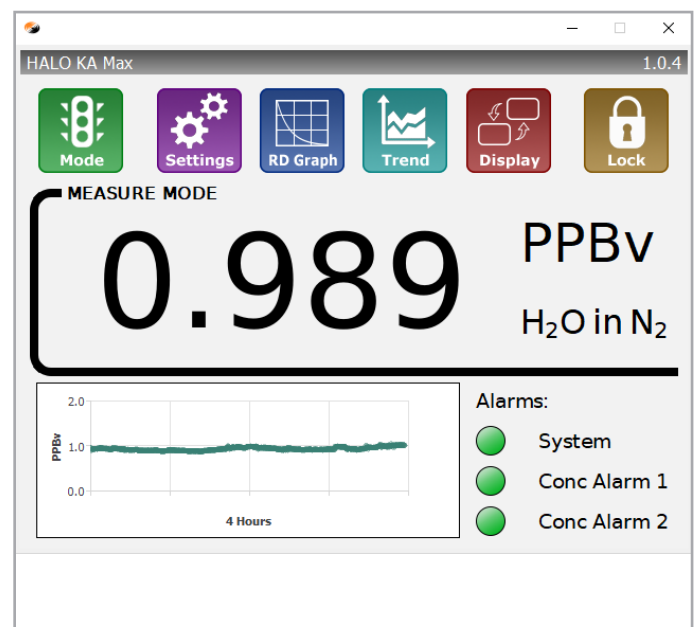


# Serani™ Max

## Analyzer Interface Software

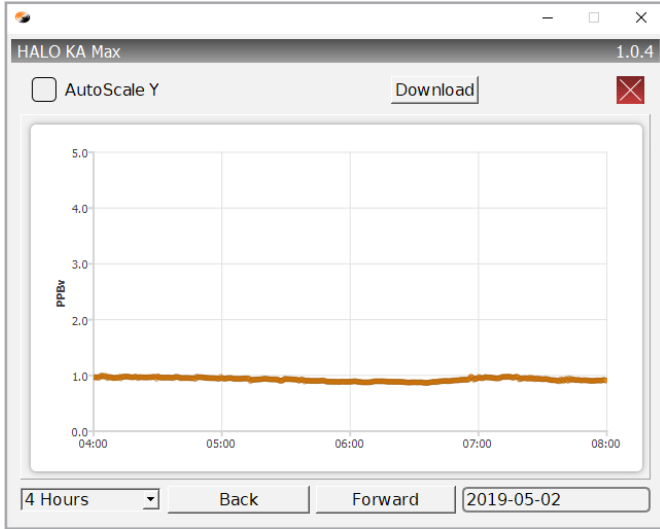

View and download the analyzer's trend data remotely from your computer or laptop.

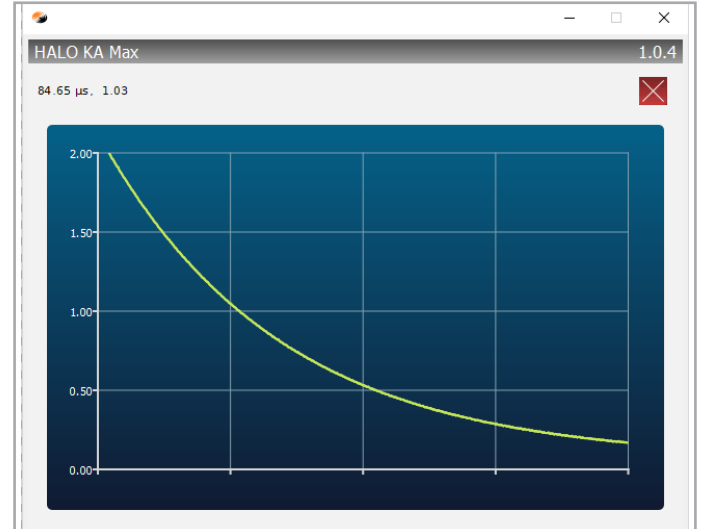

Ring-down data is streamed live to your computer to check the analyzer status remotely and in real-time.

### Available Serani™ Max Features

- ✓ Monitor Concentration in Real-Time
- ✓ View and Download Data
- ✓ Collect Verification Data
- ✓ Change Gas Type & Data Filter
- ✓ Monitor Laser and Alarm Status
- ✓ Live Streaming of Ring-Down Data

### Full List of Compatible Tiger Optics Analyzers

Analyzer Series/Family	Compatible Models
HALO	HALO Max QCL · HALO KA Max · HALO KA* · HALO 3* · HALO LP/RP/QRP* · HALO OK* · HALO H2*
Spark	Spark*
T-I Max	T-I Max · T-I Max CEM · T-I Max AIR
ALOHA H <sub>2</sub> O	ALOHA+ H <sub>2</sub> O*

\*Compatible only with analyzers that HAVE a front USB port

The Serani™ Max software is not compatible with LaserTrace and Prismatic analyzers. For HALO (3 and below, without front USB port), Tiger-i, ALOHA H<sub>2</sub>O, ALOHA+ H<sub>2</sub>O (w/o front USB port), CO-rekt and Spark (w/o front USB port) series analyzers, please see the classic Serani™ interface software.

### Distribution and Licensing

The Serani™ Max software package can be downloaded from our website, [www.tigeroptics.com](http://www.tigeroptics.com). Use of the software requires a license key, which must be purchased from Tiger Optics or an authorized representative and is tied to a specific analyzer. For more information on software licenses and pricing, or to request a demo version, please call us or contact us at [sales@tigeroptics.com](mailto:sales@tigeroptics.com).

### Tiger Optics, LLC

275 Gibraltar Road, Horsham, PA 19044  
 Phone: +1 (215) 656 4000 · Fax: +1 (215) 343 7168  
[sales@tigeroptics.com](mailto:sales@tigeroptics.com) · [www.tigeroptics.com](http://www.tigeroptics.com)

